

	TECHNICAL DATASHEET	CODICE:	N204ISO1B
		DATA:	<i>May 2016</i>
		PAGINA:	1/3

Commercial denomination

N204ISO1B

Description:

Made with EN AW-5754 UNI EN 573-3 light alloy, anodised after finishing.

Framework in rectangular tubular of 60 x 30 in galvanised hermetic anti-torsion steel.

Perimeter lower bumper in aluminium sections fitted with PVC, colour smoke grey, to prevent marking walls.

The body of the trolley is an unique piece of sheet, with vertical terminal edges flapped at the outside, without junctions and/or interstices inside, to guarantee the maximum level of hygiene and cleanness during the washing and sterilisation phase.

The top is in shaped sheet, with dripping in the anterior part, to favour the exit of the water at the exterior of the trolley.

The inside of the trolley is equipped with n.1 vertical compartment with n.10 couples of guides for the baskets; each of them can carry ISO norm baskets 600 x 400 mm.

Guides are in nylon fibre glass and stay in vertical light alloy profiles, they can be taken off with extreme facility so that they can be easily cleaned.

The trolley has only one door with closure with vertical bolt with blocking system in the upper and lower part by means of a stainless steel AISI 304 (1.4301) rod; this acts on fusions in self-lubricating nylon fibre glass and is operated by a handle.

N.1 pull handle in flat profile placed on the door.

Along the perimeter of the opening side there is a rubber gasket integrated in the structure which creates a barrier against dust (during the closing phase).

Equipped with n.4 grey anti-mark rubber wheels Ø125 mm, with supports in galvanized steel (two fix and two swivel), arranged at the corners.

Inner dim.: 465 x 625 x 1280 h mm

Outer dim.: 520 x 710 x 1500 h mm



OPTIONALS

Cod.TSA

Upper bumper in grey PVC



Cod.TB

Draw bar



Cod.CHA4

Cardholder A4




Cod.CHA5

Cardholder A5



Cod.CL

Lock with 2 keys

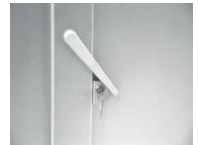

Cod.BD

Directional block on 2 wheels



Cod.RR

Pedalbrake



Cod.RRC

Centralbrake on 2 wheels

