



SHINVA MEDICAL INSTRUMENT CO., LTD

Address: Xinhua Medical Scientific Zone, Zibo New & Hi-Tech Industrial
Development Zone, Shandong Province, P.R. China

Sales Tel: +86 533 3587720 / 3587719

Email: sales@shinva.com

Website: www.shinva.com

Operating Table
STable-H8500

Electro-hydraulic Operating Table STable-H8500

STable-H8500 is a universal surgical table which offers maximum load capacity, maximum flexibility and maximum comfort. The modular design tabletop provides various configuration for OR and improves the outcomes of surgeon. With digital display on remote control, medical staff will adjust the table to desired position quickly.




Extreme Load Capacity

With a maximum weight capacity of 460kg, STable-H8500 is designed to support the lightest to the heaviest surgical patient in any surgery position, safe and durable.



Large Longitudinal slide

With longitudinal slide up to 480mm, STable-H8500 enables optimum access with the C-arm unit without reposition the patient during imaging procedure. Further, the carbon-fiber tabletop (optional) provides superior imaging quality.

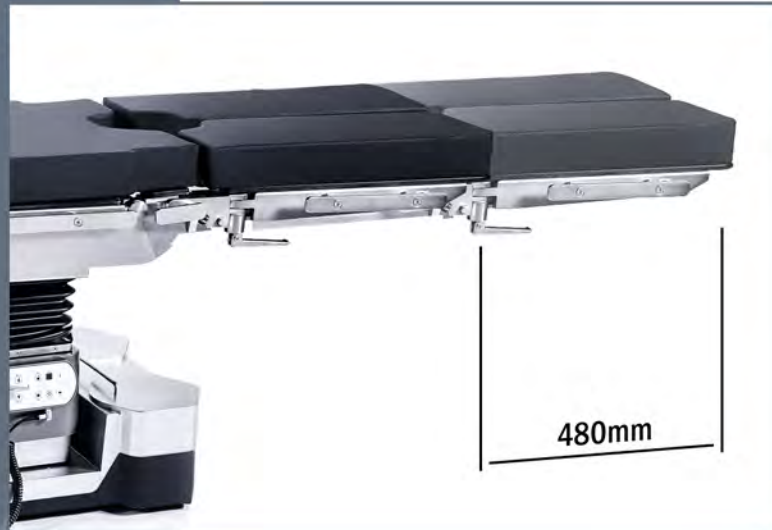

Intelligent Collision Prevention System

Built-in intelligent collision prevention system monitors the positions of the table elements during movement to avoid collision, which guarantee the safe of table.



Excellent Adjustment Design

Electro-hydraulic system provides motorized adjustment of height, Trendelenburg/Reverse trendelenburg, lateral tilt, back section, leg section individually or synchronously, flex, reflex angles, auto-leveling, floor lock and unlock, except head board.



Stable Brake With Base Bearing

Different from conventional castor bearing whole table weight design, STable-H8500 has four auto-lift castors, which means castors will extend only during table movement not all the time.



Options

Carbon fiber tabletop



Foot switch control



Orthopedic Traction Frame




Support accessories




Basic motion parameters


Tabletop elevation range
580mm-1080mm

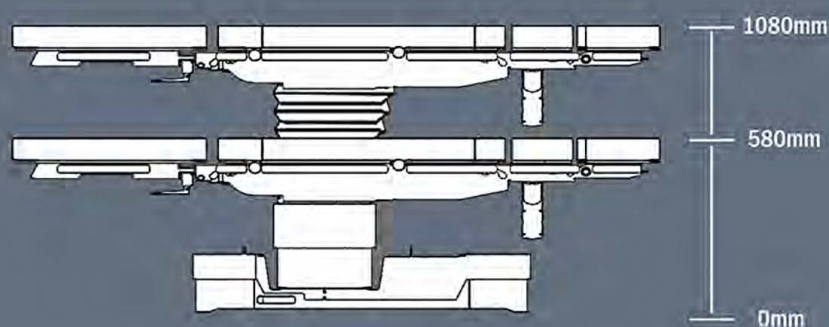


Longitudinal slide 480mm

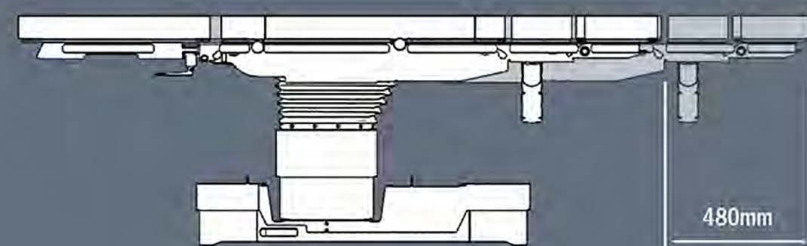


Trendelenburg 36°
Reverse trendelenburg 36°

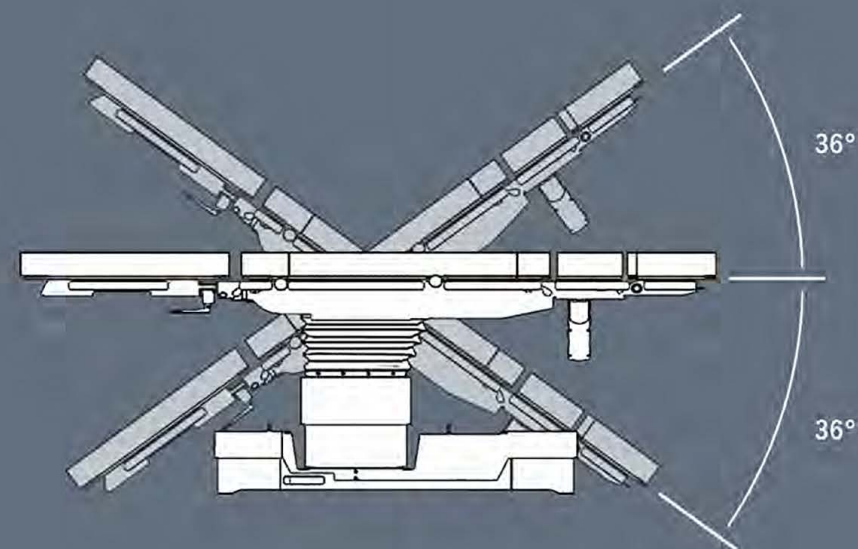


Lateral tilt left 26°
Lateral tilt right 26°

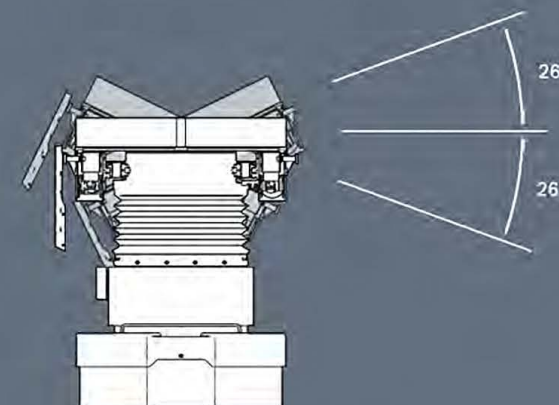


Head board up 60°
Head board down 90°

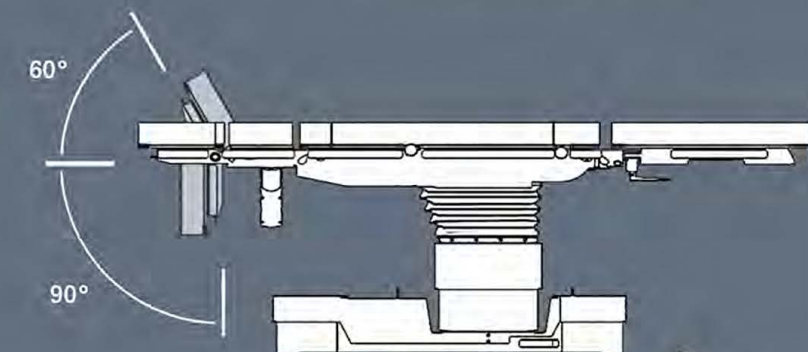

Leg plate up 80°
Leg plate down 95°

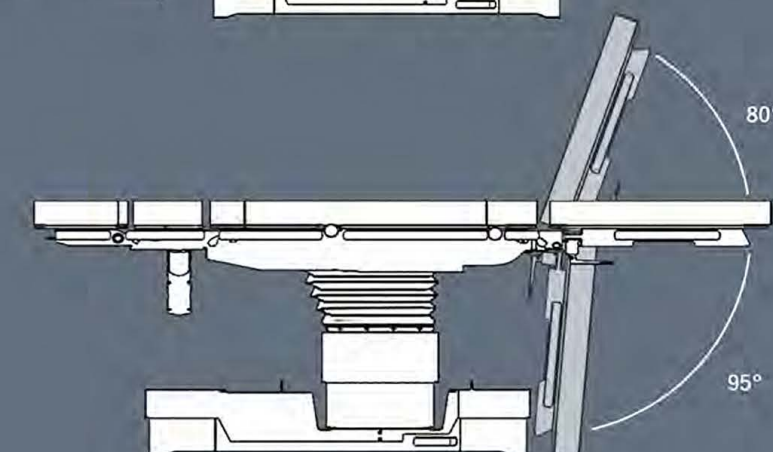

Back plate up 90°
Back plate down 45°

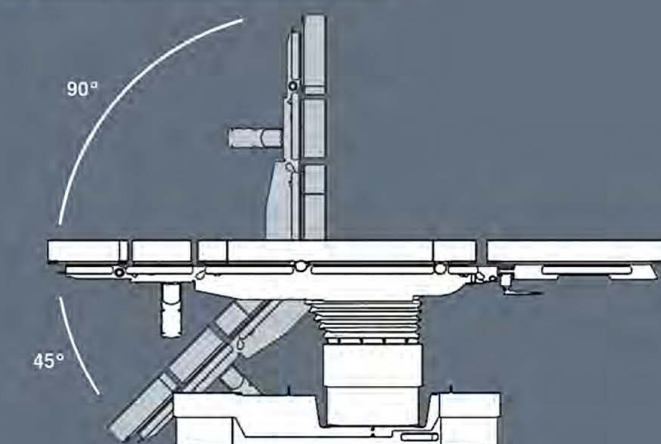

Safmed

Your total solutions partner for more than 30 YEARS

Johannesburg:
Tel: 011 201 4300
Fax: 011 507 5550

Kwa-Zulu Natal:
Tel: 082 464 3832
Fax: 031 561 5138

Western & Eastern Cape:
Tel: 021 763 3280
Fax: 021 762 2760

Free State
Tel: 082 686 8075



info@safmed.co.za



www.safmed.co.za



Safmed.SA



Safmed



Safmed South Africa