

# ADASPOR™ PLUS

HIGH-LEVEL STERILANT AND DISINFECTANT  
ISAZONE™ MOLECULE & PERACETIC ACID

## ADASPOR™ PLUS DISINFECTANT FAMILY

| DESCRIPTION                                                                                               | ACTIVITY THRU                                                                                                                                      |
|-----------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------|
| • ADASPOR™ PLUS M Disinfectant/Sterilant                                                                  | 14 Days                                                                                                                                            |
| • ADASPOR™ PLUS Concentrate Disinfectant/Sterilant                                                        | 15 Days                                                                                                                                            |
| • ADASPOR™ PLUS Pronto 5 L Disinfectant/Sterilant                                                         | 15 Days                                                                                                                                            |
| • ADASPOR™ PLUS Pronto 1 L Disinfectant/Sterilant                                                         | 15 Days                                                                                                                                            |
| • ADASPOR™ PLUS Pentadie Disinfectant/Sterilant                                                           | 5 Days                                                                                                                                             |
| • ADASPOR™ PLUS Monodie Disinfectant/Sterilant                                                            | 1 Day                                                                                                                                              |
| • ADASPOR™ PLUS Single Shot Disinfectant/Sterilant<br>ADASPOR™ PLUS Single Shot is available in 5L or 10L | For use in select<br>MEDIVATORS™ AERs                                                                                                              |
| • ADASPOR™ PLUS Test Strip                                                                                | Strips for determination of efficacy<br>and MCR (minimum recommended<br>concentration of peracetic acid in<br>ADASPOR™ PLUS Disinfectant/Sterilant |



[www.cantelmedical.eu](http://www.cantelmedical.eu)

Cantel Medical (Italy) S.r.l. - Via Laurentina, 169 - 00071 Pomezia (RM) Italia  
t: +39 06 9145399 - f: +39 06 9146099

ADASPOR™, ISAZONE™ and MEDIVATORS™ are trademarks of Cantel Medical Corp., its affiliates and related companies. All other company and product names are trademarks of their respective owners.

© 2017-2021 Cantel. All Rights Reserved. The Complete Circle of Protection diagram is a trademark owned by Cantel Medical Corp. and is subject to copyright protection owned by Cantel Medical Corp.  
LIT-00455 Rev A



# ADASPOR™ PLUS

HIGH-LEVEL STERILANT AND DISINFECTANT  
ISAZONE™ MOLECULE & PERACETIC ACID

Reprocess




[www.cantelmedical.eu](http://www.cantelmedical.eu)

# ADASPOR™ PLUS HIGH-LEVEL STERILANT AND DISINFECTANT

ISAZONE™ MOLECULE & PERACETIC ACID

Sterilization solution formulated with Peracetic Acid & Diazaadamantane (ISAZONE™ MOLECULE) with conformity to UNI EN ISO 14937 par. 5.3.1 in 10 minutes, tested under dirty and clean conditions, as required by UNI EN 14885.

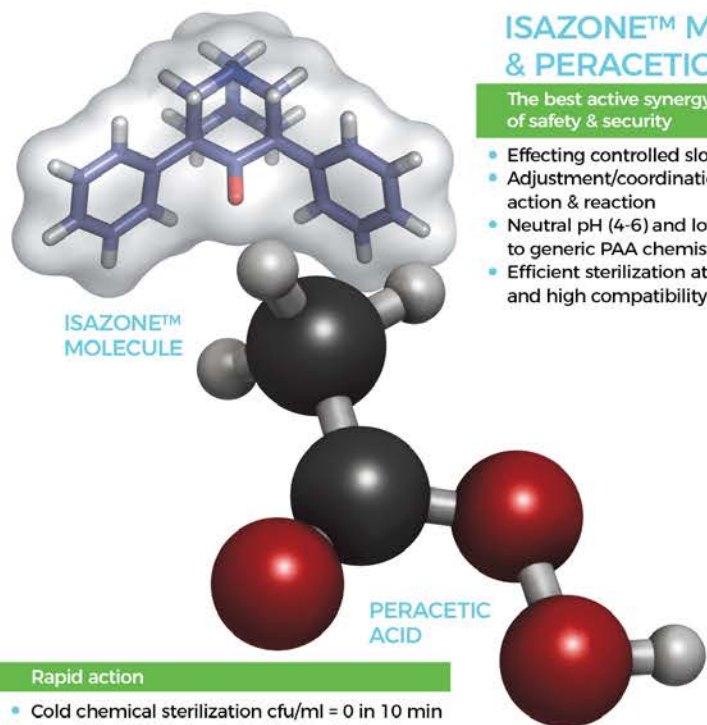

## ISAZONE™ MOLECULE & PERACETIC ACID

The best active synergy for the highest level of safety & security

- Effecting controlled slow release of active oxygen
- Adjustment/coordination of Peracetic Acid (PAA) action & reaction
- Neutral pH (4-6) and low PAA concentration in compare to generic PAA chemistries
- Efficient sterilization at ambient temperature and high compatibility



### Rapid action

- Cold chemical sterilization cfu/ml = 0 in 10 min
- Sporicidal & high-level disinfection in 5 min

### Efficiency

- Synergetic action of ISAZONE™ MOLECULE and Peracetic Acid eliminates ANY microbial form of life (cfu/ml=0), in accordance to sterilization definition (UNI EN ISO 14937- par. 5.3.1) for cold chemical sterilization
- Highest microbiological activity and widest action spectrum (sporocidal, tuberculocidal, virucidal, micobactericidal, battericidal)

## THE COMPLETE CIRCLE OF PROTECTION

As the global vanguard in infection prevention, **only Cantel delivers the Complete Circle of Protection**, a full-value, proactive partnership dedicated to helping you remove risk, streamline operational efficiencies and optimise your success.

**REPROCESS** High-level disinfection is the cornerstone of infection prevention. Reprocessing technologies from Cantel are designed to ensure patient safety by adhering to the strictest of standards for high-level disinfection, optimise workflow efficiency, and deliver versatility by supporting a wide range of scopes.



## ASSURANCE FOR EFFECTIVE REPROCESSING FOR THE PREVENTION OF INFECTIONS



### Maximal level of security

- Clinically and toxicologically tested
- Microbiological activity tests are performed in clean and dirty conditions to ensure highest safety standards for patients and operators

## COMPATIBILITY WITH THERMOLIABLE MEDICAL DEVICES

ADASPOR™ PLUS Disinfectant/Sterilant guarantees absolute compatibility with thermoliable medical devices and flexible endoscopes from different manufacturers, confirmed with millions of sterilization and high level disinfection cycles performed with ADASPOR™ PLUS Disinfectant/Sterilant worldwide. Reduction of service and maintenance costs for reusable medical devices.

### Maximum duration

- With slow controlled release due to ISAZONE™ MOLECULE
- ADASPOR™ PLUS Disinfectant/Sterilant last for 15 days after 1st usage



## ADASPOR™ PLUS TEST STRIPS

- Activity & efficacy controlled and confirmed with ADASPOR™ PLUS Test Strip
- Strips are designed for determination of MRC (Minimal recommended Concentration) 0.05% of Peracetic Acid in association with ISAZONE™ MOLECULE, ensure sterilization action of ADASPOR™ PLUS Disinfectant/Sterilant.

